



EPITOME OF PERFECTION

Accumax Bottle Top Dispenser

Accumax Bottle Top Dispenser

Where there's innovation,
there's expertise and that is
what Accumax has proved with
its marvelous Bottle Top Dispenser.

Designed after extensive research about
customer requirements and needs,
the bottle top dispenser aims to fulfill all
of your dispensing needs with unmatched
precision and convenience.




Satisfaction Guaranteed


Experience how it feels to be completely satisfied with a device which is designed and built to offer unmatched features that make using it a complete pleasure. Made from the absolute best quality raw materials, Accumax Bottle Top Dispensers promise years of durability while offering unwavering precision.



Discharge Tubes With Recirculation Valve Eliminates Sample Loss & Saves Precious Reagents & Samples

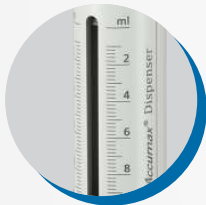

Unique Piston Mechanism Allows Cleaning And Servicing Of Piston Cylinder Without Disturbing Calibration

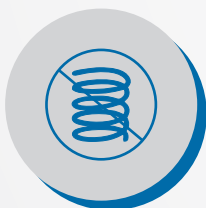

Springless Valve Design Offers High-chemical Resistance & A Longer-life With Minimal Maintenance

Volume Lock Ensures Stopping Of Any Accidental Discharge



The Device Comes With 4 Additional Adapters For Most Common Bottle Sizes For Enhanced Convenience



Versatile Telescopic Filling Tube Compatible With Different Sizes Of Reagent Bottles From 170 MM To 300 MM



A Superior Way To Work

While ensuring impeccable accuracy, Accumax's Bottle Top Dispenser is also designed to offer a superior user-experience through various features which make the device extremely easy and convenient to use. Thus, the secret behind it's wild popularity is its immense user-friendliness made possible by its highly effective features.



Discharge Tube With 360° Rotation, Makes Dispensing In Multiple Vessels Convenient



Dual Graduation Marks Enable Precise Volume Setting & Clear Visibility From Both Sides.



The Device Is Easy To Recalibrate With A Single Tool.



Calibrated As Per ISO 8655 Standard In A Lab Accredited To ISO 17025 Standard. Each Dispenser Come With Individual Calibration Report.

Ergonomic Design Enables Convenient Sampling & Low Pressure Dispensing For Even Larger Volumes



Snap-fit Mounting Of Discharge-tube For Uncomplicated Mounting



Safety Locking Feature Prevents Accidental Flow Of Reagents



Ordering Information

Accumax Bottle Top Dispenser Technical Specifications

Catalog No. With Re-circulation	Catalog No. Without Re-circulation	Volume Range ml	Increment ml	Inaccuracy (±) %	Inaccuracy (±) ml	Imprecision (±) %	Imprecision (±) ml
ADR - 2.5	AD - 2.5	0.25 to 2.5	0.05	0.6	0.015	0.2	0.005
ADR - 05	AD - 05	0.5 to 5	0.1	0.5	0.025	0.2	0.01
ADR - 10	AD - 10	1 to 10	0.2	0.5	0.05	0.2	0.02
ADR - 25	AD - 25	2.5 to 25	0.5	0.5	0.125	0.2	0.05
ADR - 30	AD - 30	2.5 to 30	0.5	0.5	0.15	0.2	0.06
ADR - 50	AD - 50	5 to 50	1.0	0.5	0.25	0.2	0.1
ADR - 60	AD - 60	5 to 60	1.0	0.5	0.3	0.2	0.12

Accessories

Description	Cat#	Color
Bottle Cap 28mm	ADA-01	
Bottle Cap 38mm	ADA-02	
Bottle Cap 40mm	ADA-03	
Bottle Cap 45mm	ADA-04	
Calibration Tool	ADA-05	
Recirculation Tube	ADA-06	
Telescopic Fill Tube Asm 2.5-60ml	ADA-07	



Accumax Lab Devices Pvt. Ltd.
www.accumaximum.com



ISO 9001:2015
ISO 13485:2016

* OEM request can be entertained. Specifications can be changed without notice for quality improvement.